



Invacare® | **PEGASUS METRO**

The perfect mix of comfort,
safety and reliability



Whether you enjoy a quick trip to meet friends for food or taking a short journey to the local shops, the **Pegasus METRO** is the ideal choice to take you there.

The new, ergonomically designed **Pegasus METRO** has been designed from the ground up to focus on key aspects important to both new and seasoned scooter users. As a result, it is equipped with numerous features to ensure users experience a safe, easy and enjoyable ride.

The **Pegasus METRO** is the perfect mix of comfort, safety and reliability!



Controlled and sturdy drive

The suspension system on the compact **Pegasus METRO** has been designed to navigate a variety of everyday surfaces and terrains, allowing you to relax and enjoy a controlled, sturdy drive in comfort. It's powerful motor, combined with 11" pneumatic tyres, will take you smoothly and safely over obstacles or uneven ground.

Safety as standard

In addition to a unique lighting system, the **Pegasus METRO** is equipped with an electronic speed reduction feature that automatically reduces speed round bends, ensuring safe, comfortable manoeuvres at all times. Furthermore, a two-step disengaging lever prevents the scooter from free-wheeling if it is knocked out of drive mode. These features, combined with extensive German independent TÜV laboratory testing, ensure that safety comes as standard with the *Invacare* **Pegasus METRO**.

Simplicity and serviceability

By using an array of quality components, we have made the **Pegasus METRO** both reliable and easy to service. Simply remove the top cover to access batteries, electronics, motor and cables. Electronics and the motor are protected from corrosion and water which safeguard the vital workings of the scooter against dirt preventing any unnecessary malfunctions.



Features and options




Finite tiller movement

Users can easily adjust the tiller with a lever to suit their needs.



Ergonomic steering

Our ergonomically designed steering system makes for easier, more responsive control.

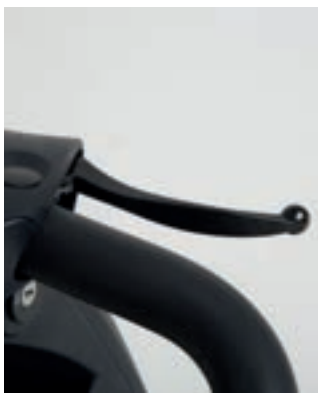

Hand brake

Ensures immediate braking if required.



Brake light

Clear, bright brake lighting warns people that the scooter is slowing down.



Disengaging lever

Two step disengaging lever prevents the scooter from accidentally free-wheeling.














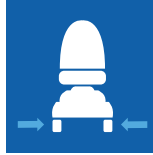


Auto-speed reduction




Automatic speed reduction reduces speed around bends for added safety.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

	Backrest width	Seat depth	Seat to floor height	Backrest height	Driving unit width	Total length	Empty weight
<i>Pegasus Metro</i>	 510 mm	 Comfort and Deluxe: 470 mm Premium: 460 mm	 Min. 440 mm Max. 510 mm	 Comfort: 500 mm Deluxe: 490 mm Premium: 630 mm	 660 mm	 1270 mm	 110 kg


	The heaviest part	Max. user weight	Battery capacity	Turning diameter	Max. climbable obstacle height	Ground clearance	Max. safe slope
<i>Pegasus Metro</i>	 57 kg	 136 kg	 2 x 12 V/50 Ah (C20) AGM	 260 0mm	 80 mm	 78 mm	 6°

	Max. driving distance	Max. driving speed	Motor capacity
<i>Pegasus Metro</i>	 54 km	 6 / 10 / 12km/h	 1. Peak performance 6 km/h: 240 W / 500 W ¹ 10 km/h: 240 W / 600 W ¹ 12 km/h: 550 W / 1500 W ¹

Colour of shrouds

Marble White	Sapphire Blue	Jade Green	Luxor Beige	Crystal Copper
				
Carmine Red	Iconic Chrome	Diamond Silver	Electric Blue	Jasper Red
				

Upholstery

Black Vinyl


Invacare Australia Pty Ltd:
ABN 45 074 676 378
PO Box 5002, 1 Lenton Place, North
Rocks, NSW 2151, Australia.
Freephone: 1800 460 460,
Fax: 02 8839 5353
Email: orders@invacare.com.au
Web: www.invacare.com.au

Invacare New Zealand:
PO Box 62-124, 4 Westfield
Place, Mt. Wellington,
Auckland, New Zealand.
Freephone 0800 468 222,
Freefax 0800 807 788
Email: sales@invacare.co.nz
Web: www.invacare.co.nz

Invacare Asia Ltd:
PO Box 5002, 1 Lenton Place,
North Rocks, NSW 2151,
Australia.
Phone: +61 2 8839 5330,
Fax: +61 2 8839 5331
Email: asiasales@invacare.com
Web: www.invacare.com

Trusted Solutions,
Passionate People

