



Invacare® | **PEGASUS PRO**

Comfort, safety and reliability
with added performance



Ideal for longer trips into town or afternoons driving around the park with your grandchildren, the new ergonomically designed **Pegasus ^{PRO}** is a powerful upgrade to our standard **Pegasus ^{METRO}** model.

With enhanced battery performance, superior suspension and large 12" wheels, the new **Pegasus ^{PRO}** is ideal for users who need that extra performance.

The **Pegasus ^{PRO}** blends this extra performance with all of our existing comfort, safety and reliability features to make it the ideal choice for users who demand more.



Advanced suspension, smoother drive

The advanced suspension system on the **Pegasus ^{PRO}** has been designed to bring better performance and give a smoother ride over a variety of surfaces and terrains. In addition, the **Pegasus ^{PRO}** features a powerful motor and 12" pneumatic tyres to ensure that whether you are in the city or the country, you get where you need to go in comfort and style.

Performance with security and safety

The **Pegasus ^{PRO}** has been subjected to German independent TÜV laboratory testing, ensuring that users can be confident that they are driving a scooter that adheres to the latest European Medical Device Standards.

In addition, an on-board electronic speed reduction system automatically reduces speed round bends, ensuring users are confidently making safe, comfortable manoeuvres. Furthermore, a two-step disengaging lever prevents the scooter from free-wheeling if it is knocked out of drive mode.

Quick and easy serviceability

The **Pegasus ^{PRO}** has been built using an array of quality components, which have enabled us to make it both reliable and easy to service. The electronics and motor have additional protection to guard against water and corrosion, safeguarding the vital workings of the scooter against any unnecessary malfunctions. The top cover can also be easily removed to check the motor, batteries, electronics and cables.



Features and options



Finite tiller movement

Users can easily adjust the tiller with a lever to suit their needs.



Ergonomic steering

Our ergonomically designed steering system makes for easier, more responsive control.



Hand brake

Ensures immediate braking if required.



Brake light

Clear, bright brake lighting warns people that the scooter is slowing down.



Disengaging lever

Two step disengaging lever prevents the scooter from accidentally free-wheeling.



Auto-speed reduction

Automatic speed reduction reduces speed around bends for added safety.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

	Backrest width	Seat depth	Seat to floor height	Backrest height	Driving unit width	Total length	Empty weight
<i>Pegasus Pro</i>	 510 mm	 Comfort and Deluxe: 470 mm Premium: 460 mm	 Min. 440 mm Max. 510 mm	 Comfort: 500 mm Deluxe: 490 mm Premium: 630 mm	 650 mm	 1320 mm	 136 kg

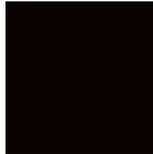
	The heaviest part	Max. user weight	Battery capacity	Turning diameter	Max. climbable obstacle height	Ground clearance	Max. safe slope
<i>Pegasus Pro</i>	 approx: 62 kg	 160 kg	 2 x 12 V/75 Ah (C20) AGM	 2750 mm	 100 mm	 100 mm	 10 °

	Max. driving distance	Max. driving speed	Motor capacity
<i>Pegasus Pro</i>	 52 km	 6 / 10 / 12 / 12.8 / 15 km/h	 1. Peak performance 6 km/h: 240 W / 500 W ¹ 10 km/h: 240 W / 600 W ¹ 12 km/h: 550 W / 1500 W ¹ 12.8 km/h: 550 W / 1600 W ¹ 15 km/h: 550 W / 1600 W ¹

Colour of shrouds

Marble White	Sapphire Blue	Jade Green	Luxor Beige	Crystal Copper
				
Carmine Red	Iconic Chrome	Diamond Silver	Electric Blue	Jasper Red
				

Upholstery

Black Vinyl


Invacare Australia Pty Ltd:
ABN 45 074 676 378
PO Box 5002, 1 Lenton Place, North
Rocks, NSW 2151, Australia.
Freephone: 1800 460 460,
Fax: 02 8839 5353
Email: orders@invacare.com.au
Web: www.invacare.com.au

Invacare New Zealand:
PO Box 62-124, 4 Westfield
Place, Mt. Wellington,
Auckland, New Zealand.
Freephone 0800 468 222,
Freefax 0800 807 788
Email: sales@invacare.co.nz
Web: www.invacare.co.nz

Invacare Asia Ltd:
PO Box 5002, 1 Lenton Place,
North Rocks, NSW 2151,
Australia.
Phone: +61 2 8839 5330,
Fax: +61 2 8839 5331
Email: asiasales@invacare.com
Web: www.invacare.com

Trusted Solutions,
Passionate People

